

OPTIFOAM GENTLE WITH A LIQUITRAP CORE Super Absorbent with a Gentle Touch





Highly breathable film acts as a waterproof secure barrier to bacteria

Foam layer absorbs exudate to help minimize pooling of exudate at the wound surface

Super absorbent layer absorbs and retains exudate in the dressing minimizing the risk of maceration


Absorbent net facilitates absorption of exudate into the super absorbent layer and helps prevent leakage

Silicone adhesive minimizes trauma to the wound and reduces pain upon wear or removal



Item No.	Description	Pkg.
MSC2333EP	3 x 3"	10/bx, 100/cs
MSC2344EP	4 x 4"	10/bx, 100/cs
MSC2366EP	6 x 6"	10/bx, 100/cs
MSC2377EP	Sacrum, 7 x 7"	5/bx, 40/cs
MSC2399EP	Sacrum, 9 x 9"	5/bx, 25/cs

Borderless - cut to size

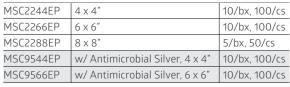

BENEFITS

- » Super absorbent Liquitrap core draws in large amounts of heavy exudate, converts it into a gel and locks it away
- » Low profile design increases flexibility and comfort
- » Five unique layers provide added protection

FOLLOW US fin BLOG



Medline Industries, Inc. One Medline Place, Mundelein, IL 60060

Medline United States 1-800-MEDLINE (633-5463) medline.com | info@medline.com

Medline Canada 1-800-396-6996 medline.ca | canada@medline.com

Medline México 01-800-831-0898 medlinemexico.com | mexico@medline.com